


LABCONCO LOW PROFILE FUME HOOD SYSTEM

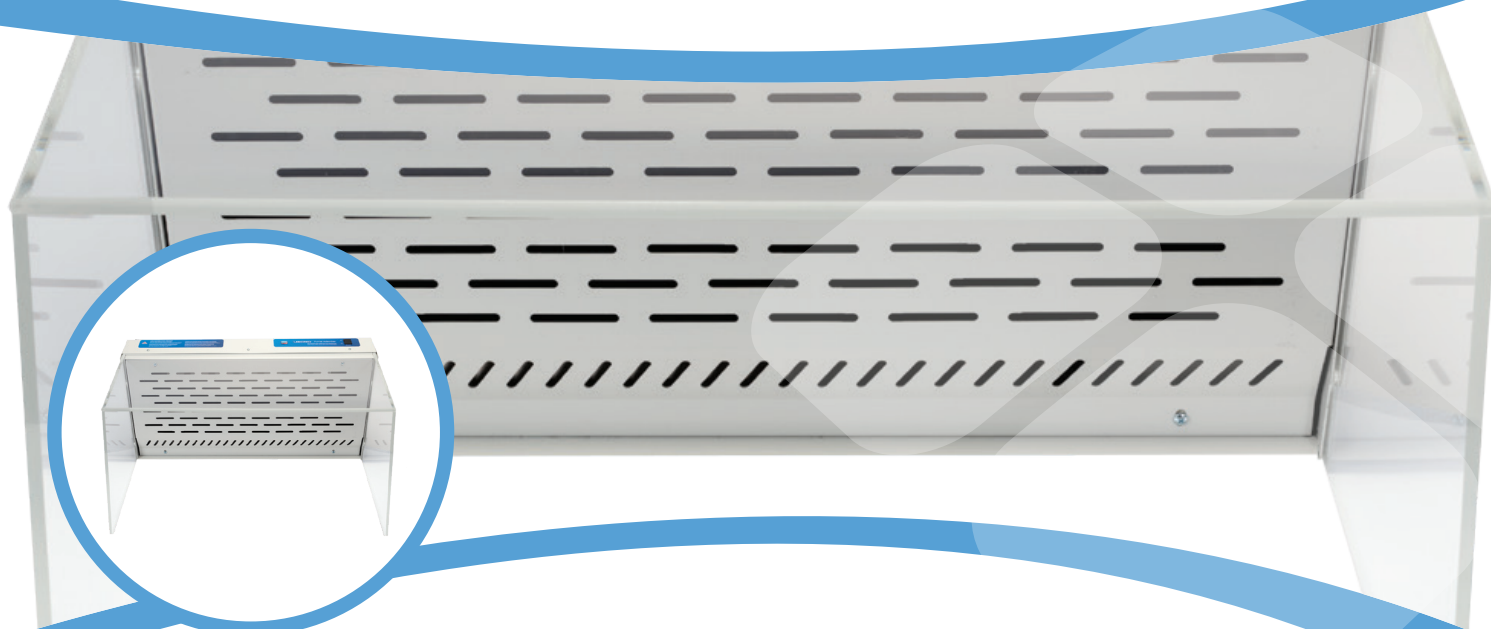

Labconco Low Profile Fume Hood Systems are geared for light duty histology coverslipping and staining areas. Activated carbon media adsorb organic odors and noxious vapors while clean air is recirculated into the laboratory. These enclosures are low profile, completely portable and require no duct work. They easily accommodate shallow counter tops or those with overhanging cabinets. Each unit includes either an activated carbon, formaldehyde or ammonia filter.

Features

- Built-in motorized impeller
- Adsorber has a white, dry powder epoxy-coated steel housing with perforated baffle for uniform airflow
- Canopy is a clear, one-piece 3/8"-thick, shatter-resistant acrylic

Benefits

- Carbon filtration traps vapors so that safe, filtered air is returned to the laboratory
- No outside duct work is required
- These portable, low-profile enclosures fit easily on countertops with overhanging cabinets

Product Information

SY5014-A Labconco Low Profile Fume Hood System

Related Accessories


SY5019 Activated Charcoal Filter

Call us today.

800-783-9424
info@avantik-us.com
avantik-us.com



LABCONCO LOW PROFILE FUME HOOD SYSTEM

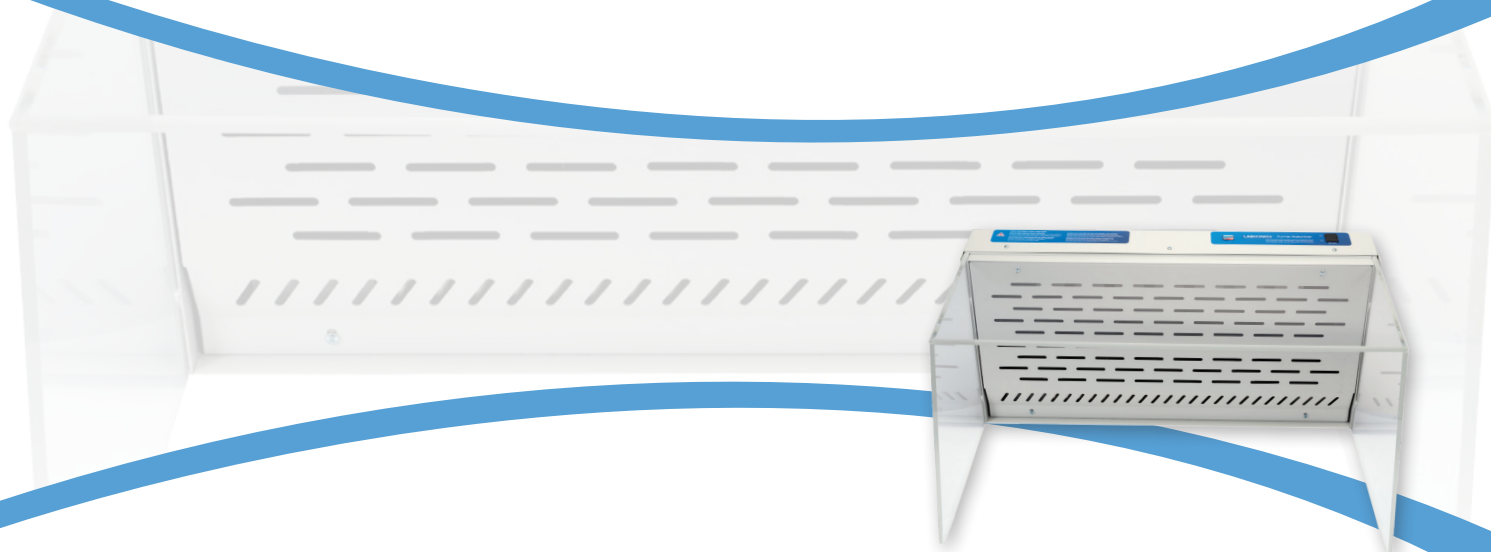

Perfect for the coverslipping and staining area of any histology lab. The Labconco Low Profile Fume Hood System uses charcoal media to absorb annoying odors and noxious vapors such as xylene and toluene. When operated as directed, levels of these toxic vapors fall below OSHA standards. The Labconco Low Profile Fume Hood System will easily fit on any countertop, even with overhanging cabinets. Lightweight and portable, the Fume Hood System requires no outside venting.

SY5014-A Labconco Low Profile Fume Hood System

Height	15.4 in.
Width	30 in.
Depth	23.4 in.
Weight	62 Lbs.

Technical Specifications

Power	115V AC, 50/60 Hz, 1 amp AC (230V AC optional for export)
Blower	100 CFM (free air) Backward Curved Impeller
Noise Level	Weighted Noise Level 68 dB at 1 Meter

Standard Fume Inclusions

- Shatter resistant enclosure (1/4" thick one piece, clear acrylic)
- Two replaceable filters (AC activated carbon filters)
- Wiring
- 8 foot, grounded cord & plug

Available Options

- 230 volt operation - 50/60 Hz, 1 amp AC (Requires an electrical plug, not included.)
- Replacement charcoal filters (6 pounds of type AC activated charcoal media 26"W x 10.5"H)

Call us today.

800-783-9424
info@avantik-us.com
avantik-us.com

